

kuraray

CLEARFIL™ REPAIR Multi Purpose

A System for Repairing All Intra-oral Restorations

A Porcelain Repair System That Works – And Lasts

- Simple, fast, reliable repair of chipped or fractured porcelain, composite resin, or metal oxide restoration
- Very high bond strength even to metal oxides, like zirconia
- No toxic Hydrofluoric Acid needed
- Long, stable shelf-life (two years from manufacturing)
- Re-glazing porcelain and composite restorations



KURARAY AMERICA, Inc.

600 Lexington Avenue, New York, NY 10022 | 800.879.1676 | kuraraydental.com | ©2008 Kuraray America, Inc.

Clinical Applications

- Repairing fixed prosthetic restorations fabricated from porcelain, metal oxides (including zirconia), hybrid ceramic, or composite resin (e.g., inlays/onlays, crowns, bridges, and the superstructure of implants)
- Re-glazing fixed intra-oral restorations (e.g., cured direct/indirect composite resin, provisionals, acrylic appliances, glass ionomers)


Repairing fixed prosthetic restorations fabricated from porcelain, ceramics, hybrid ceramics or composite resin.

<p>1 Roughen the surface with a diamond point, rinse and dry.</p> 	<p>2 Apply K-etchant GEL on the surface, leave in place for 5 seconds, rinse and dry.</p> <p>Zirconia and alumina prosthetic restorations do not require this step.</p> 	<p>3 If the surface contains metal, apply ALLOY PRIMER.</p> 
<p>4 Apply CLEARFIL™ CERAMIC PRIMER to the surface and dry.</p> 	<p>5 If the surface contains metal, apply CLEARFIL™ ST OPAQUER and light-cure.</p> 	<p>6 Place and light-cure the selected composite resin (e.g. CLEARFIL MAJESTY™ Esthetic).</p> 



CLEARFIL™ CERAMIC PRIMER

- Does not require toxic Hydrofluoric Acid (HF)
- Use with any porcelain, composite, hybrid ceramics, and metal oxide ceramics (alumina and zirconia)
- Fast and easy application with one component (No mixing needed)

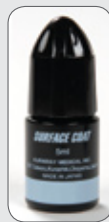


SURFACE COAT™

- Colorless, transparent gloss that “feels good” to your patient
- Glazed surface resists staining
- High gloss that can protect the restoration surface



ALLOY PRIMER

- Increases bond strength to gold and semi-precious metals
- Simple to use, eliminates tin-plating
- Intra- and extra- oral use

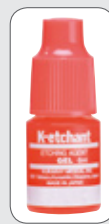

K-etchant GEL

- 40% phosphoric acid etching gel
- For etching teeth, porcelain, and cured composite for adhesive restorations

ITEM#	PRODUCT NAME	CONTENTS
#1973KA	CLEARFIL™ REPAIR Multi Purpose	K-etchant GEL 1 bottle, 6 ml CLEARFIL™ CERAMIC PRIMER 1 bottle, 4 ml ALLOY PRIMER 1 bottle, 5 ml CLEARFIL™ ST OPAQUER (US) 1 syringe, 4.0 g SURFACE COAT™ 1 bottle, 5 ml Accessories 2 Brush tip handle (black and white) Disposable brush tips 50 Normal type / 20 Flat type 1 Mixing dish 1 Light blocking plate 1 Paper pad

REFILL		
#013KA	K-etchant GEL	1 bottle, 6 ml
#2550KA	CLEARFIL™ CERAMIC PRIMER	1 bottle, 4 ml
#064KA	ALLOY PRIMER	1 bottle, 5 ml
#636KA	CLEARFIL™ ST OPAQUER L	1 syringe, 4.0 g
#637KA	CLEARFIL™ ST OPAQUER US	1 syringe, 4.0 g
#263KA	SURFACE COAT™	1 bottle, 5 ml

ACCESSORIES		
#902	Mixing dish	1
#965	Brush tip handle (White)	1
#966	Brush tip handle (Black)	1
#968	Disposable brush tip	1 (50 tips)
#1010KA	Disposable brush tip (flat type)	1 (20 tips)

Warranty: Kuraray Medical Inc. will replace any product that is proved to be defective. Neither Kuraray Medical, Inc. nor the distributor accepts liability for any loss or damage, direct, consequential or special, arising out of the application or use of or the inability to use these products. Before using, the user should determine the suitability of the products for the intended use and the user assumes all risk and liability whatsoever in connection therein. CLEARFIL™ REPAIR Multi Purpose and Kuraray products are available only through Kuraray Authorized dealers. CLEARFIL MAJESTY™ and CLEARFIL™ are trademarks of Kuraray Co., LTD. © 2008 Kuraray America, Inc.